

Cat. No. FYD009-1ML

100 bp YEA Ladder DNA Marker II

Ready-to-use DNA ladder

Recommended storage condition: -20°C

Loading volume:	5~10 µl/well
Size range:	100-3000 bp
Concentration:	60 ng/µl
Package size:	60 µg/1 mL

100 bp DNA Ladder Plus 1500 bp, 2000 bp and 3000 bp, ready-to-use

The Yeastern Biotech's 100 bp YEA Ladder DNA Marker II (FYD009-1ML) contains 13 discrete fragments ranging from 100 bp to 3000 bp. It is supplied in a ready-to-use format containing blue tracking dyes, and mixed at varied concentrations.

For best results, load 10 µL of the marker per well. It allows sizing and concentration estimation of DNA fragments on agarose gels generated by PCR or restriction digestion.



www.yeastern.com


yeastern@yeastern.com

t / 886-2-26953922 f / 886-2-26953979

Cat. No. FYD009-1ML

100 bp YEA Ladder DNA Marker II

Ready-to-use DNA ladder

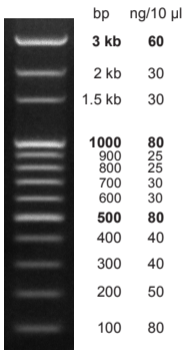

Loading volume: 10 μ l/well
1.5%, 1 \times TAE agarose gel



Yeastern Biotech Co., Ltd.